

SML™ Dorsal

SML™ Dorsal utilizes precision-milled upper and lower occlusal splints to open the vertical relationship while holding the mandible in the desired forward position. Interlocking buccal inclines allow mandibular lateral movement, alleviating muscle trismus. The inclines, while keeping the mandible from dropping open during sleep, enable the patient to open and close normally while wearing the appliance. The SML™ Dorsal is ideal for patients with existing crown and/or bridge work.



SML™ Dorsal	\$499	16 days in lab
-------------	-------	----------------

TAP® - Thornton Adjustable Positioner

The **TAP®** series of products are mandibular advancement devices for the treatment of snoring and sleep apnea. TAP products fit comfortably in the mouth, provide more room for the tongue, and allow the lips to close. Interchangeable hooks accommodate Class I, II, and III bites and increase vertical distance between the upper and lower teeth.

TAP 3®	\$485	16 days in lab
dreamTAP®	\$565	16 days in lab



Myerson EMA® - Elastic Mandibular Advancement

The **EMA®** is a customized removable appliance for the non-invasive treatment of snoring and Obstructive Sleep Apnea. It is designed to advance the mandible and open the bite to allow for less restricted airflow during sleep.

Myerson EMA®	\$400	16 days in lab
--------------	-------	----------------




Herbst

The **Herbst** is Medicare-approved as a sleep apnea and orthodontic appliance. The one-piece telescopic mechanism advances the mandible into a forward position allowing for an improved airway. Write "medicare patient" on your prescription if applicable.

The Snap Fit Herbst uses Adams clasps, while the Sleep Herbst uses ball clasps. The Snap Fit Herbst is designed for short teeth with low retention.

Sleep Herbst	\$456	16 days in lab
Snap Fit Herbst	\$821	16 days in lab

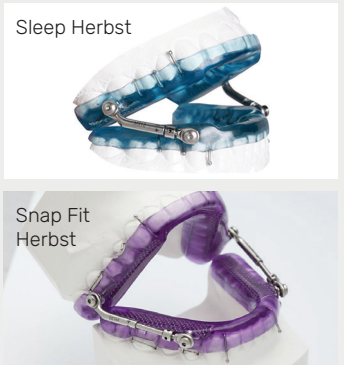

Order Today www.roedentallab.com/Removables